
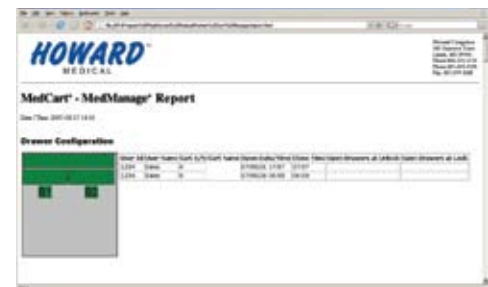
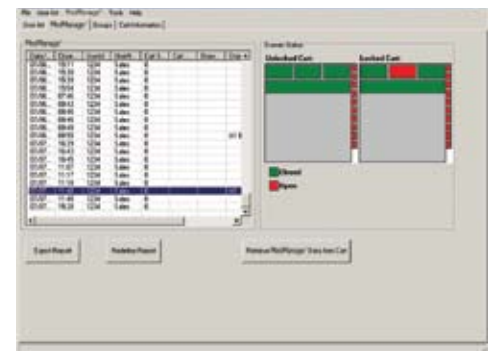


CART SOFTWARE PACKAGE

1. MedDisplay software allows the cart's user interface to be carried out via the notebook or a touchscreen on the cart itself. This permits you to:
 - Open/close the drawers electronically
 - Provide access to the narc lock electronically
 - Indicate individual drawers that are open
 - Display how much time until the drawers will lock electronically
2. MedManager software provides communication between the medical cart and a personal computer. The MedManager software can be used to maintain a centralized database of all user-IDs and the MedAudit records, as well as enable cart programming. This includes:
 - Add/delete user IDs (up to 200 IDs)
 - Change timeout on automatic locking
 - Change the cart's name
 - Run MedAudit (audit trail)
 - Search for carts on the network
 - Add temporary user IDs
3. The MedAudit software maintains a record of all transactions that occur when the cart is locked and unlocked. This audit trail creates an event record that comprises the following information:
 - Date and time the cart is unlocked
 - Date and time the cart is re-locked
 - User ID that unlocked the cart
 - Storage drawers and bins that are opened while the cart is unlocked; and storage drawers and bins that remained open after cart re-lockedThe MedAudit stores up to 500 event records using a first-in-first-out recording method. Thus, the most recent events will always be in the software. The report can be reviewed on the monitor; as an option, the report can be output to a personal computer for archiving. When the data is transmitted, the cart's ID will be sent with the report so that the data can be appropriately correlated to that cart.

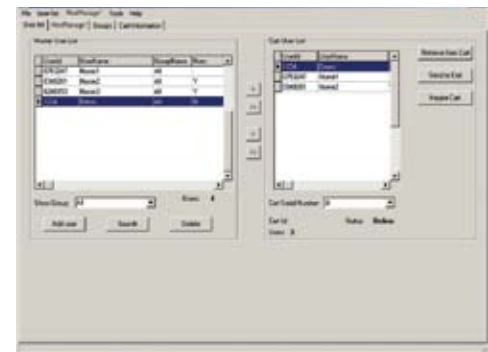

HOWARDTM
MEDICAL

Visit us online at Howard-Medical.com or call us at 877.856.6441